


THE COMPLETE

qPCR Workflow Checklist

A complete workflow reference for planning your assay setup before you order

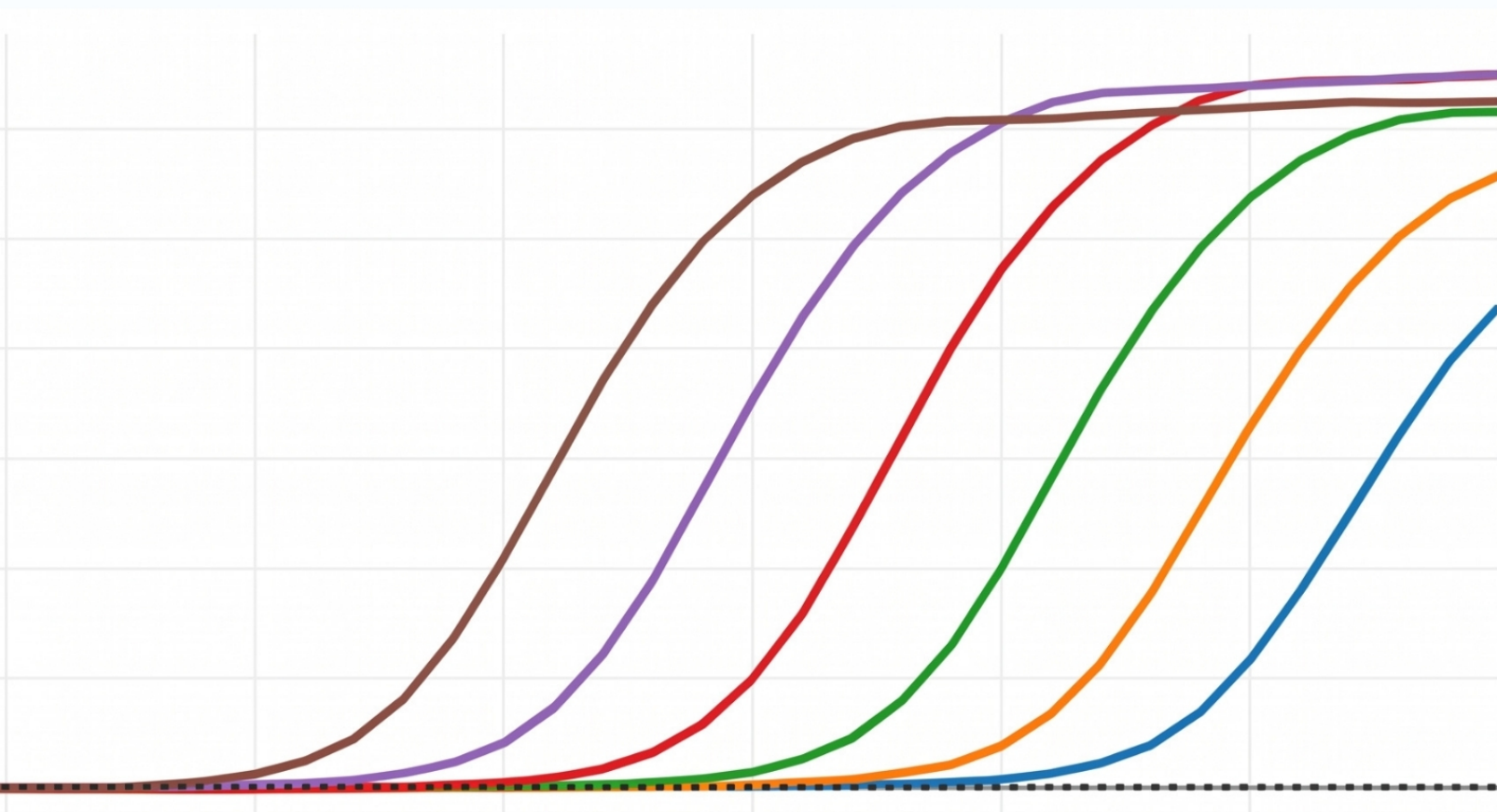

The Lean Lab Playbook is a practical operational guide for lab managers, directors, and senior scientists. Nothing herein constitutes billing, legal, or compliance advice. Consult your billing team and legal counsel regarding reimbursement policies and payer contracts specific to your lab. BioPathogenix Products are RUO for Research Use Only. Not for use in diagnostic procedures.

WHAT YOU NEED

A pre-order checklist for a complete qPCR workflow. Work down it before you order — the gaps are usually in the last two groups.

COMPONENT

WHERE BIOPATHOGENIX COVERS IT



01 SPECIMEN COLLECTION & TRANSPORT

- | | |
|--|-----------------------------------|
| ☐ Collection devices | Specimen collection supplies |
| ☐ Transport medium for your matrix | Amies, Cary-Blair and other media |
| ☐ Specimen bags and sharps containers | Consumables catalogue |
| ☐ Laboratory PPE | Gloves, coats, masks, shoe covers |



02 NUCLEIC ACID RELEASE

- | | |
|---|--|
| ☐ Extraction chemistry | DNA / RNA extraction reagents |
| ☐ Automation, if throughput needs it | NucleoMini (BPX-16E) · BPX Extraction Instrument |
| ☐ Or an extraction-free route | DNA Direct-Detect · Duo Direct-Detect |
| ☐ A way to demonstrate the stage ran | qCHECK™ external positive control |



03 REACTION SETUP

- | | |
|---|---|
| ☐ Master mix | Master Mix Developed In-House |
| ☐ Primer and probe set | Optimised sets, supplied with the kit |
| ☐ Reaction vessel, plate and consumables | Lyophilized 96- and 384-well, breakaway |
| ☐ A written protocol | Kit user manual and protocols |



04 AMPLIFICATION & DETECTION

- | | |
|--|---|
| ☐ Real-time PCR instrument | nucleomini extraction machine |
| ☐ Enough channels for the multiplex level | OctaQuant (coming soon) provides eight optical channels |
| ☐ A dye set your platform can read | qPLEX™ panels in matched configurations |
| ☐ A cycling program | Supplied with the kit protocol |



05 CONTROLS & VERIFICATION

- | | |
|---|-----------------------------|
| ☐ External positive control through extraction | BPX™qCHECK™ Quality Control |
| ☐ Internal control | Carried in-panel |
| ☐ Instrument verification | RNase P Verification Plate |
| ☐ Characterised material for method work | Biobank · contrived samples |

Not sure which route fits? Turn over.

order@biopathogenix.com · support@biopathogenix.com · (859) 444-5660

For Research Use Only. Not for use in diagnostic procedures.

**No order minimums
No contracts**

THREE ROUTES

What changes if you run it differently.

Groups 01, 04 and 05 stay the same whichever route you take. Groups 02 and 03 are where the configuration choice lands.

Standard extraction

Best fit: Established workflows handling mixed sample types

- 02** Extraction reagents plus a purification instrument if throughput calls for it
- 03** Master mix, primer and probe set, plate and consumables — prepared and distributed

Direct-Detect

Best fit: Hands-on time and throughput pressure

- 02** Two reagents. No purification chemistry, no purification instrument
- 03** Unchanged — master mix, primer and probe set, plate and consumables

Lyophilized plates

Best fit: Cold storage limits, field and remote sites

- 02** Unchanged — extraction chemistry and instrument as normal
- 03** The plate is the setup. No master mix to prepare, nothing to distribute

COMMONLY MISSED

-  **A control that reaches extraction**
An amplification-only control passes while extraction underperforms. The negative looks clean.

-  **A channel for the internal control**
A four-target panel needs five channels, not four. The control is easy to forget when counting.

-  **Instrument verification**
Nothing else in the workflow reports on whether the block is cycling and reading uniformly.

-  **Matrix-matched material**
Method work needs material in the matrix the assay will actually see, not a convenient substitute.